



# Versto

Lightweight Adaptable  
**POWERED WHEELCHAIR**

**VERMEIREN**  
we care for you

# Verso

## Lightweight Adaptable **POWERED WHEELCHAIR**

Verso is versatile and agile powered wheelchair. Its lightweight design makes it a perfect wheelchair for users on the move who would like to extend their mobility without compromising on comfort.

Thanks to our extensive expertise in options, Verso offers its user a remarkable adaptability to its user. Crafted by merging the durability of a manual wheelchair frame with the advanced features of two motorized rear wheels, Verso is a perfect choice for a comfortable, lightweight powered wheelchair. The lithium battery can be removed for easy charging to enhance convenience in your daily use.

**VERMEIREN**  
we care for you







# STANDARD PACKAGE

Removable and adjustable armrest

Height adjustable push handles – 5 positions

Height adjustable backrest with 4 straps for tension adjustment

Breathable 3D NET cushion upholstery

Fold away and removable footrest, adjustable in length

Removable lithium battery, equipped with safety lock



Adjustable seat depth from 42 to 48 cm

# EASY TO USE

Simple, yet perfect. Verso has it all. Dismount your wheelchair in a few seconds. All parts are easy to carry and not too heavy for easy transport in a car.

Assembly is easy with our Click&Drive system. No tools needed, just connect your battery and you're ready to drive! Thanks to our removable battery, you can charge your chair whenever and wherever you like.



**Verso is compatible with aircraft transport!**

Just switch off the battery and your wheelchair can fly. Always refer to the user manual for transport instructions.

# EASY TO TRAVEL

A unique feature for Verso is the detachable motor and battery block. Detach your wheelchair from its battery and wheel block to allow for easy storage and transport. Both the frame and the “powered part” of Verso have an additional wheel to make sure no parts can be damaged by detaching. They also allow for easy manoeuvring.

The seat can easily be folded like you would any of our other top-selling manual wheelchairs. Perfect for transport in the booth of the car or storage in your home.



# LIGHTWEIGHT

**31,3 kg**

Total weight



**9,9 kg**

Transport weight  
(without cushion, footrest and armrest)



**16,2 kg**

Heaviest part





**We care about you... and your back.** Verso is a very light powered wheelchair with only 16,2 kg for the motorised part and 9,9 kg transport weight for the main frame. Easy to lift without any strain.





# DRIVING EXPERIENCE



## Controller

The P&G controller is just how we like it at **VERMEIREN**, no-nonsense. With a battery status indicator, speed adaptor and a comfortable joystick, this controller is easy to use. Basic programming is also possible to personalise your driving experience.



## Up to 21 km

Experience an extended driving range with the 17.5 Ah battery. Select the battery that suits your longing for an escape.



## 6 km/h

With 6 km/h Verso allows you to drive at the perfect speed for comfortable conversation and walks with your friends and family.



# BATTERY - MOTORS

The Lithium-ion battery (available in 12,8 Ah or 17,5 Ah) has a battery status indicator and is lightweight. Easily unlock your battery with your own anti-theft key to charge your battery anywhere you like. No time to charge? An additional battery is available as an option.

The battery is conveniently located at the back of your wheelchair for easy access when you need to remove it for charging or to quickly check your battery status.

With 2 x 250 W motors and 12" PU black tires, Verso is competing with others in its category. A great and durable choice for a powered wheelchair for everyday use.

Removable  
Li-ion  
battery

12,8 Ah - 24 V  
17,5 Ah - 24 V

2 x 250 W  
motors





CAUTION \* Store in cool and dry place! \* Do not expose the product to fire, water or moisture

100% 80% 50% 20%

A close-up photograph of a bicycle handlebar. The handlebar is black and features a bright blue rectangular plate with the word 'Verso' in a stylized, embossed font. The plate is secured with two black screws. In the background, the black frame of the bicycle and a portion of a black tire with a tread pattern are visible, slightly out of focus. The lighting is dramatic, highlighting the textures and colors of the components.

# DESIGN & COLOUR

With a sleek design and fun colour accents, Verso can be adapted to your preferences whilst having a modern and attractive look. The colour accents serve a double purpose as they allow you to easily locate the levers of the powered part. For additional safety we included some reflectors.

# COLOURS

Colours are a personal experience, just like your Verso. Choose the colour combination that suits your style the best.

Select your frame colour between glossy white and glossy black and play with colour accents to make it your style.

## Available colour accents





# HIGHLY CONFIGURABLE

Adapt your Verso to you. Verso is a highly configurable chair which combines the large and extensive option range from **VERMEIREN** and our expertise in powered wheelchairs. Choose whatever suits you best!



# TWO in ONE

CONVERT **Verso** INTO A MANUAL WHEELCHAIR

## KIT 24" wheels


Verso offers a quite unique option for people who are in need of both a manual and powered wheelchair. With our 24" wheel kit you can easily change between manual or powered driving.

**Plug and drive!**



# OPTIONS

With a very large range of options, Verso is perfect to suit individual needs as well as the needs of the carer. Configure the chair you need with our wide selection of options.

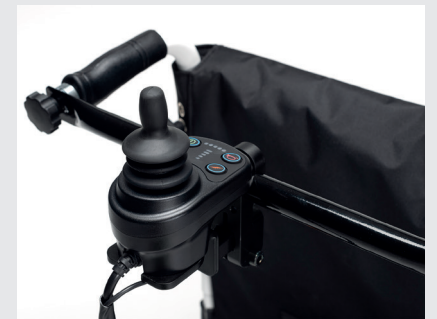

*For more options please refer to your local order form.*



# TECHNICAL DATA

■ Total width	seat width + 17cm
■ Total length	1078 mm
■ Seat height	521 mm
■ Seat width	440 - 460 - 480 - 500 mm
■ Seat depth	420 to 480 mm
■ Backrest height	400 - 420 - 440 - 460 - 480 mm

■ Batteries	12,8 Ah - 24 V   17,5 Ah - 24 V
■ Motors	2 x 250 W
■ Speed*	6 km/h
■ Driving range*	15 km   21 km
■ Weight	31,3 kg
■ Maximum user weight	130 kg

\*Measurements in mm and kg. Deviations in sizes and weight (15 mm and 1,5 kg) are possible  
Subject to technical alterations. Images for illustration purposes only. Certain features may be optional.



# Versto

Lightweight Adaptable  
**POWERED WHEELCHAIR**

## **VERMEIREN Group**

Vermeirenplein 1-15  
B-2920 Kalmthout  
Belgium

[www.vermeiren.com](http://www.vermeiren.com)  
[info@vermeiren.com](mailto:info@vermeiren.com)



For more information, contact your local dealer