




TriGO^o



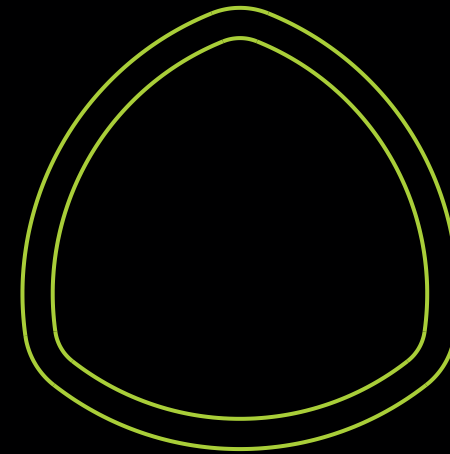

Why triangular tubes?

Our research and development teams have perfected the triangular tubes that form the basis of your **Trigo**. Apart from giving your wheelchair a modern and unique look, they have the advantage of making your Trigo lighter while maintaining optimal strength.

TRIGO S



Our wheelchair is available in two versions. First, there is **Trigo S**, an active and foldable chair with integrated footrests into the frame. The design of Trigo S is active and sporty. This configuration makes the wheelchair even lighter, while maintaining maximum strength.





TRIGO T



There is also a second model, **Trigo T**. This wheelchair has been designed with the same modern look. The foldable frame has retractable and removable footrests. Its design makes this model ideal for easy transport.

DESIGN

VERMEIREN has designed **Trigo** as an ultra-modern wheelchair to suit everyone's taste. Therefore, your Trigo can be fully personalised upon ordering; from the frame to the wheels, including the armrests and a range of more than 30 possible colours. In short, there are almost as many different Trigo wheelchairs as there are users!





INFINITELY ADJUSTABLE

Apart from the configurations to be chosen upon ordering (size, type of backrest, armrests, frame, options,...), your **Trigo** remains a dynamic chair. A large number of adjustments can be made on each chair, such as wheel camber, seat height,...



ADJUSTABLE WHEEL CAMBER FOR EACH WHEELCHAIR

Thanks to an ingenious system, the wheel camber of your **Trigo** can be altered directly on your chair! The camber can be altered in the desired position of 0° , 2° , 4° or 6° (depending on configuration). Your Trigo will be delivered with the camber as chosen upon ordering.





POSITIONING OF THE WHEELS

To adapt **Trigo** to the user and to make it more or less 'active', the position of the rear wheels can be adjusted directly on the frame of your wheelchair. The rear wheels can be adjusted forward or backward into 8 different positions. Your Trigo is delivered with its wheels in the position as asked upon ordering.



BACKREST INCLINATION

There are 3 possible backrests for **Trigo**: a fixed backrest, an inclinable backrest, and an inclinable and foldable backrest. For the second and third options, the backrest inclination can be easily adjusted on the chair frame.



AVAILABLE IN CHILDREN'S SIZES

The **Trigo** wheelchair has been designed in such a way that it is suitable for adults as well as children. For this reason, Trigo is available in seat widths 30 cm up to 50 cm. Even for the smaller sizes, a wide variety of accessories and options is available.



OPTIONS

PUSH HANDLES



PH 01 Fixed push handles (12 cm)
PH 02 Short, fixed push handles (8 cm)



PH 03 Fixed, foldable push handles



PH 04 Height-adjustable push handles
(inside backrest tubes)



PH 05 Height-adjustable push handles
(behind backrest tubes)



PH 06 Foldable and height-adjustable
push handles (behind backrest tubes)



Push handles **PH 03** and **PH 06**
are fully foldable.

OPTIONS

ARMRESTS AND CLOTHING PROTECTORS



**B82-S or
B82-C FIXED**
Fixed, adjustable
clothing protectors,
aluminium or carbon



B82 FIXED
Fixed, adjustable clothing
protectors, aluminium with
mudguard



Bi5
Height-adjustable
armrests



TUBE
Armrests,
retractable to the rear
(only with fixed backrest)



B82SL-A
Long, removable
armrest, aluminium



B82L-A
Long, removable and
adjustable clothing protectors
with mudguard, aluminium



B82S-C
Removable and
adjustable clothing
protectors, carbon



B82-C
Long, removable and
adjustable clothing protectors
with mudguard, carbon

Available in angles of 70°, 80° and 90°



Monobloc tube in aluminium



Monobloc, footplate in aluminium, foldable



Monobloc, footplate in carbon, foldable



Monobloc, with higher footplate



Monobloc, with retractable footplate



Separate footplates

OPTIONS

FOOTRESTS



70°



80°



90°



BA6 footrests, shortened by 5cm



OPTIONS

WHEELS

Front wheels



Standard black
3", 4", 5", 6"



Transparent
3", 4", 5", 6"



Aluminium silver
3", 4", 5"



Aluminium blue
3", 4", 5"



Aluminium red
3", 4", 5"



Carbon
4", 5"



LED
4", 5"

Rear wheels



Standard wheels 18"
18"



Design wielen
20", 22", 24", 26"



Spirit Wheels Aluminium
24"



Spirit wheels red
24"



Spirit wheels blue
24"



Off-road wheels
24"

Measuring unit: inches



Right Run black tyres
available with Design
and Spirit wheels



Marathon Plus tyres
available with Design
and Spirit wheels

OPTIONS

BRAKES



Push or pull brakes



Sport brakes



Traditional brakes with retractable handles (B40)

OPTIONS

OPTIONS AND ACCESSORIES



**Crutch holder
left or right**
B31L or B31R



**Anti-tipping wheel
left or right**
B78L or B78R



Travel wheels
B86



**Stabilisation bar
for backrest**
B90



Neoflex belts
Consult all possible
solutions on our website
www.vermeiren.com



Pouch at the rear



**Multi-adjustable
headrest**
B47



**Side supports left
or right (small size)**
L04S L or L04S R



Drum brake
grey wheels
24" x 1/3 8"



Spoke protectors
Consult all possible
models on our website
www.vermeiren.com

More options and accessories to be found on our order forms.

MULTIPLE COLOUR OPTIONS!



C4M - Bright Red Mat
C4G - Bright Red Gloss



C5M - Black Mat
C5G - Black Gloss



C8M - Yellow Mat
C8G - Yellow Gloss



C29M - Bright Grey mat
C29G - Bright Grey Gloss



C53M - Orange Mat
C53G - Orange Gloss



C55M - Inox Mat
C55G - Inox Gloss




C56M - Blue Grey Mat
C56G - Blue Grey Gloss



C59M - Champagne Mat
C59G - Champagne Gloss



**C71G - Metallic
Deep Blue**

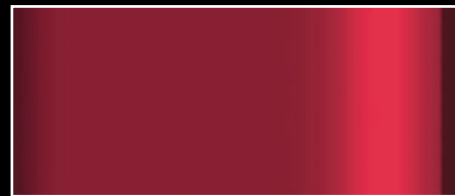

**C72G - Metallic
Cherry Red**



**C73G - Metallic
Bright Green**



C76M - Metallic White Mat
C76G - Metallic White Gloss



**C78G - Soft
Blue Gloss**



C79M - Metallic Bright Blue Mat
C79G - Metallic Bright Blue Gloss



C85M - Grass Green Mat
C85G - Grass Green Gloss



C87M - Melon Orange Mat
C87G - Melon Orange Gloss



C88M - Gold Mat
C88G - Gold Gloss
(with surcharge)



C89M - Copper Mat
C89G - Copper Gloss
(with surcharge)



**Dozens of other
colours at your request.**
(with surcharge)

TECHNICAL DATA



seat width : **300 up to 500 mm** (steps of 20 mm)

total width (camber 0°) : **460 up to 660 mm**

seat depth : **320 up to 500 mm** (steps of 20 mm)

seat height : **380 up to 510 mm** (according to configuration)

total weight : **starting from 12 kg** (according to configuration)

maximum user weight : **70 or 130 kg** (according to configuration)

possibility of foldable backrest : **yes**

backrest height : **250 up to 510 mm** (according to configuration)

front wheels : **3", 4", 5", 6"**

rear wheels : **18", 20", 22", 24", 26"**

wheel camber : **0° up to 6°** (directly adjustable on the chair)

Available upholstery: **Nylon, Nylon PLUS, 3D Net**

Dimensions in mm, weight in kg. Measurement tolerance: 15 mm, 1,5 kg, 1,5°



VERMEIREN

N.V. VERMEIREN N.V.

Vermeirenplein 1-15
2920 Kalmthout
Belgium

Tel : +32 (0)3 620 20 20
Fax : +32 (0)3 666 48 94
@ : info@vermeiren.com

www.vermeiren.com

All rights reserved, including translation. • No part of this manual may be reproduced in any form what so ever (print, photocopy, microfilm or any other process) without written permission of the publisher, or processed, duplicated or distributed by using electronic systems.

Edition: September 2020.

Images shown only for illustration purposes.
Certain features may be optional. Subject to technical alterations.

Visit our website the for the most up to date information: www.vermeiren.com



TriGO^o

